


Poly-L-Lysine Coating (PLL)



The NEXTERION® poly-L-lysine surface enhances cell attachment and adhesion of biomolecules to our glass surface using ionic interactions. This coating is based on long chains of positively charged L-lysine polymers. Produced in ISO class 5 clean room conditions using a standardized process following a stringent quality control system, NEXTERION® poly-L-lysine coated substrates are available in both standard and custom formats.

Product Information

Coating Chemistry

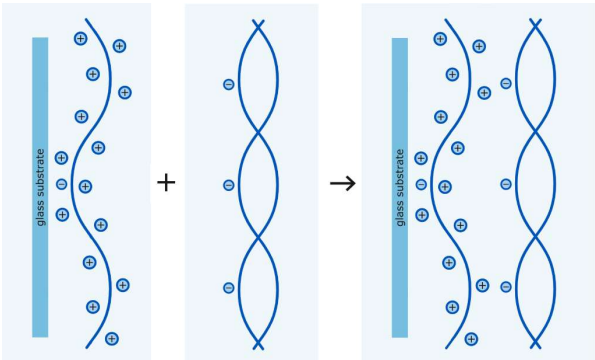

NEXTERION® poly-l-lysine coating chemistry

Shelf Life

Six months for sealed packages at room temperature.


Immobilization Method

Ionic interactions

Probe Types

- BACs, PACs, YACs
- Oligonucleotides ≥ 40 mers
- cDNA
- PCR products



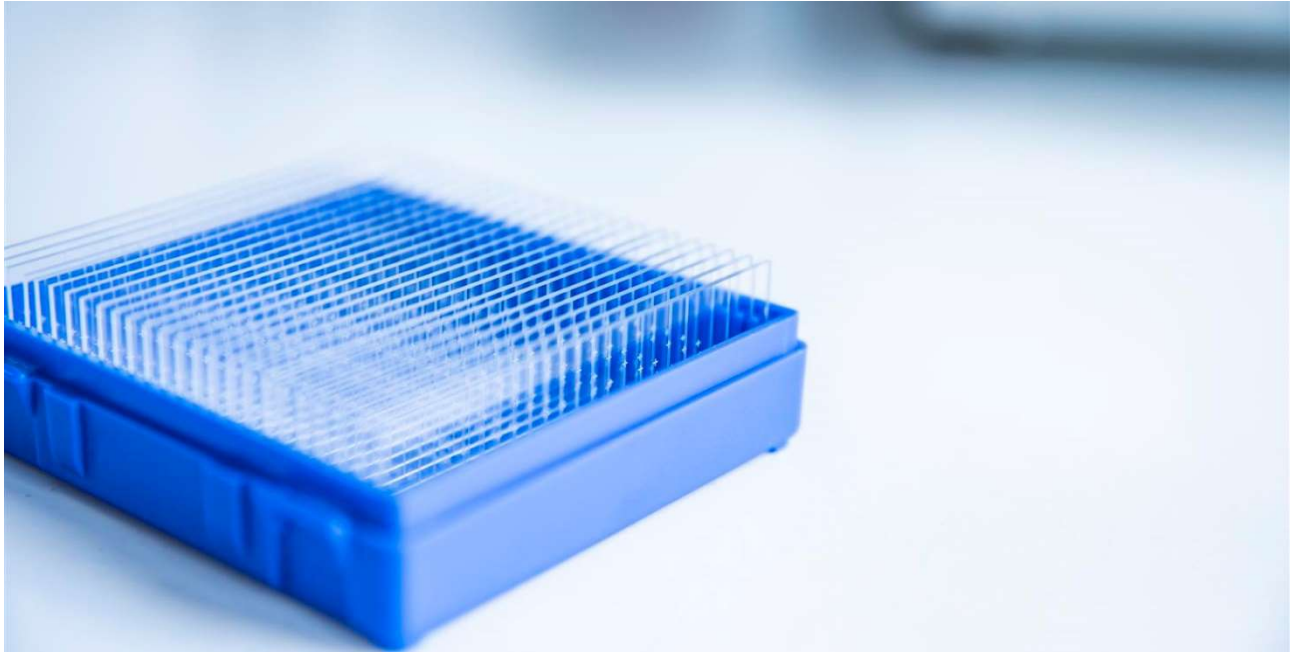

Advantages

Material

- High-quality borosilicate glass
- Alternative substrate materials can be offered

Formats

- Standard sizes (slide format)
- Customized dimensions and thicknesses

Structuring

- Pre-scoring
- Hydrophobic coating for multiplexing

Markings

- Barcodes e.g. Code128, DataMatrix-Code
- Logos
- Position markings and fiducials

Quality

- Proprietary thin-film deposition process optimized by SCHOTT
- Excellent intra- and inter-lot reproducibility
- Physical and functional quality control
- ISO class 5 clean room production
- Relevant processes in place for diagnostic company needs

Supply Forms

Product	Size (mm)	Thickness (mm)	Pieces per pack
Slide PLL	75.6 x 25.0	1.0	25
Customized PLL	Variable	0.1 - 2.5	Variable

SCHOTT Technical Glass
Solutions GmbH
Otto-Schott-Straße 13
D-07745 Jena, Germany
Phone +49 (0)3641/681-4066
Fax +49 (0)3641/681-4970
info.nexterion@schott.com

www.schott.com/nexterion

SCHOTT
glass made of ideas